G3FERRARI®

BY TREVIDEA

MOD.:

G30021

MANUALE D'USO

USER MANUAL

IT

EN

FR

Technical model: EB522

Bilancia pesabimbi - Baby scale

CRESCOBENE

www.g3ferrari.it

ITA: INDICE

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ENGLISH: SAFETY INSTRUCTIONS

The following are important notes on the installation, use and maintenance. Save this instruction manual for future reference. Use the equipment only as specified in this guide. Any other use is considered improper and dangerous. Therefore, the manufacturer cannot be held responsible in the event of damages caused by improper, incorrect or unreasonable use.

Before use, ensure the equipment is undamaged. If in doubt, do not attempt to use it and contact the authorized service center. Do not leave packing materials (i.e. plastic bags, polystyrene foam, nails, staples, etc.) within the reach of children as they are potential sources of danger; always remember that they must be separately collected.

In case of failure and / or malfunction do not tamper with the unit. To repair the appliance, always contact the service center authorized by the manufacturer and request the use of original spare parts. Failure to comply with the above may compromise the unit safety and invalidate the warranty.

- Do not handle or touch the appliance with wet hands or bare feet.
- Do not expose the appliance to harmful weather conditions such as rain, moisture, frost, etc.. Store the unit in a dry place.
- Children should be supervised to ensure that they do not play with the appliance; do not allow incapable persons to use the product without supervision.
- Place it on flat and stable surface to ensure a correct operation. Keep at a distance of 2 meters from microwave ovens in operation, which can cause malfunctions due to electromagnetic interference.
- Keep the unit out of the reach of children and away from heat sources and steam.
- After finishing weight operation, empty the scale.
- The buttons should be pushed gently and only with fingers.
- If the display does not light up, it could be a battery problem: check that they are correctly installed and they are not exhausted.
- Not suitable for commercial use. Household use only
- Not immerse in water or other liquids and avoid splashing of water; clean with damp cloth only.
- Do not subject the product to strong impacts, serious damage may result.
- Do not use as a shelf for supporting objects.
- **WARNING: this is not a diagnostic instrument, but is suitable to follow the weight variation of the baby. For more accurate and professional weighting, refer to your pediatrician.**

G30021

INSERT BATTERIES

The compartment for the battery is located under the scale.

The scale works with 1 x CR2032 lithium battery.

- Open the door by pulling the lock tab.
- Remove any battery by gently pulling on the metal tab on the edge of your battery.
- Insert the battery with the plus sign (+) facing up.
- To insert, push the battery into its slot until it stops.
- Close the cover.

CAUTION

- Always use lithium batteries CR2032 type.
- Be sure to insert the battery according to the correct polarity as otherwise they may leak and damage the scale. Such damage is not covered by warranty.
- When you are going not to use the unit for a long time, it is advisable to remove the battery to prevent the run-out.
- Do not heat, disassemble or throw in fire.
- Batteries are collected separately. Dispose of exhausted batteries correctly, as required by national regulation.

OPERATING THE UNIT

- **WARNING: this is not a diagnostic instrument, but is suitable to follow the weight variation of the baby. For more accurate and professional weighting, refer to your pediatrician.**

Before use, select the correct measuring unit by using the selector under the scale.

Weighing

- Place the scale on stable and flat surface.

- Press  and wait until on display appears “0.000 kg”

- Place the baby on the scale with care and attention.

- The scale is equipped with a HOLD system:

1) Even if the baby slightly moves, the automatic hold function can fix the weight.

2) If the baby moves too much, the weight cannot be visualized. In this case you can press the “HOLD” key, and after 3 seconds the weight will be displayed.

- The automatic weighting (option 1) is better because it is more accurate. Therefore try to calm down the baby without touching him, in case he moves too much.

- After fixing the weight, it will be displayed for 20 seconds, then the scale will automatically switch off.

Tare function

If you don't want to place the baby directly on the scale, but on a towel, you can use the tare function.

- Press  and wait until on display appears “0.000 kg”

- Place the towel on scale.

- After the word HOLD on display stops flashing, press again  and wait until on display appears “0.000 kg”

- Weigh the baby.

Note: if you turn on the scale AFTER having placed the towel it is not necessary to tare.

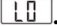
Select the measurement unit

Using the key in the base, the unit of measurement can be changed from kilograms (kg) to pounds (lb) to ounces (oz). The display shows "kg" or "lb" "oz", depending on the unit chosen.

Exceeding load

The capacity of the scale is up to 20 Kg. Once exceeding that value, the display shows “EEEE”.

Low battery indication

When the battery is running low, the display shows .

Replace the battery with a new one, following the steps shown in the "Insert batteries" section.

SWITCH-OFF The switch-off is automatic: once the weighing value is displayed, after 20 seconds the scale is switched off automatically.

CLEANING Clean the unit with a damp cloth and soft and non-abrasive cleaners only.

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WASTE OF ELECTRIC AND ELECTRONIC EQUIPMENT

The product is made of non-biodegradable and potentially polluting substances if not properly disposed of; other parts can be recycled. It's our duty to contribute to the ecological health of the environment following the correct procedures for disposal. The crossed out wheeled bin symbol indicates the product complies with the requirements of the new directives introduced to protect the environment (2002/96/EC, 2003/108/EC, 2002/95/EC, 2012/19/EC) and must be properly disposed at the end of its lifetime. If you need further information, ask for the dedicated waste disposal areas at your place of residence. Who does not dispose of the product as specified in this section shall be liable according to the law.

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To find the nearest service center, please call the telephone number 0541 694246, fax number 0541 756430 or contact us at assistentatecnica@trevidea.it asking for the service centers supervisor.

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TREVIDEA S.r.l. si scusa inoltre per eventuali errori di stampa.
TREVIDEA S.r.l.